



Reju™

Instructions for Use





Treating with Silk'n Reju for the First Time

The skin should be, clean, dry and free of any powder, gel or cream.

- Remove Silk'n Reju *APPLICATOR* from its cradle in the *BASE UNIT*.
- Press the *POWER ON/OFF SWITCH*. The green *READY/TREATMENT MODE I - INDICATOR LIGHT (I)* will start blinking, and shortly afterwards will be lit steadily.
- Your Reju is ready for treatment.
- Place the *TREATMENT SURFACE* on your skin. The *LED* light will automatically activate.
- Apply light pressure while slowly moving the Reju *APPLICATOR* in small circular motions. The *TREATMENT SURFACE* temperature will slowly increase.
- Perform the treatment for 5-7 minutes on each area.
- It is recommended to apply a small amount of a moisturizing cream on the treated skin areas after each session and to massage it gently in circular motions.

Table of Contents

Before You Start	2
What is Silk'n Reju?.....	2
Silk'n Reju device description.....	2
Package contents.....	3
Safety with Silk'n Reju	3
When to avoid using Silk'n Reju?.....	3
Reducing the risk of injury.....	4
Possible side effects using Reju.....	4
Protect your Silk'n Reju.....	4
Skin Treatment the Reju Way	5
The benefits of Silk'n Reju.....	5
Silk'n Reju facial treatment areas.....	5
Plan your Silk'n Reju treatment for best results	
First Steps with Silk'n Reju	5
Silk'n Reju setup.....	5
Silk'n Reju battery charging.....	5
Treating with Silk'n Reju for the first time.....	6
What to expect when treating with Silk'n Reju?.....	7
After treating with Silk'n Reju.....	7
Maintenance of Silk'n Reju	8
Cleaning Silk'n Reju.....	8
Storage.....	8
To protect the environment.....	8
Removing the rechargeable batteries.....	8
Troubleshooting	9
Customer Service	10
Specifications	11
Labeling	12

Before You Start

Before using Silk'n Reju for the first time, please read this Manual in its entirety. Pay particular attention to sections on device use procedures, device operation, and after-use procedures.

What is Silk'n Reju?

Silk'n Reju is an FDA cleared hand held light emitting diode (LED) device for facial treatments of periorbital wrinkles, reducing pore size and temporary improvement of the appearance of the skin. The light emitted by Silk'n Reju stimulates deep skin tissue and achieves facial rejuvenation.


Typical areas for Silk'n Reju treatment are the neck, the cheeks, the forehead, and around the eyes.

Silk'n Reju device description

Silk'n Reju device consists of a *BASE UNIT*, an *APPLICATOR* and an *ADAPTOR*.

Silk'n Reju *APPLICATOR* is a hand held unit that contains *TREATMENT SURFACE* at the *APPLICATOR* tip, *ON/OFF SWITCH* and *INDICATOR PANEL*.

Silk'n Reju *BASE UNIT* is used to charge the *APPLICATOR* batteries.

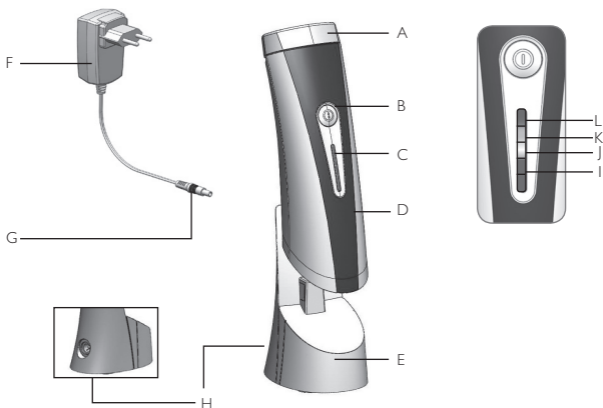

A	Treatment surface	G	Plug
B	On/Off Switch	H	BASE UNIT socket
C	Indicator Panel	I	Ready/Treatment – ON/OFF Green indicator light
D	APPLICATOR	J	System Error - Yellow indicator light
E	BASE UNIT	K	Battery charging - Orange indicator light
F	Adaptor	L	Battery full – Green indicator light

Package contents

Upon opening the Silk'n Reju package, you will find the following:

- Reju BASE UNIT and APPLICATOR
- Adaptor
- Post treatment cream
- User manual

Safety with Silk'n Reju

To protect your skin, a temperature sensor is embedded in the treatment tip of the Silk'n Reju. This sensor ensures the skin's surface will not exceed 41°C / 106°F.

When to avoid using Silk'n Reju?

Certain conditions may temporarily prevent the use of Silk'n Reju. **DO NOT USE Silk'n Reju** if any of the following conditions apply to you:

- If you have dark brown or black spots, such as large freckles, birth marks, moles or warts on the area to be treated.
- If you have eczema, psoriasis, lesions, open wounds or active infections, such as cold sores in the area to be treated. Wait for the effected area to heal before using Silk'n Reju.
- If you have abnormal skin conditions caused by diabetes or other systemic or metabolic diseases.
- If you are currently or have recently been treated with Alpha-Hydroxi Acids (AHAs), Beta-Hydroxi Acids (BHAs), Retin-A®, topical retinoids or azelaic acid.
- If you have been treated with Accutane® (isotretinoin) within the past 6 months.
- If you have been on a steroid regimen within the past 3 months.
- If you have a history of herpes outbreaks in the area of treatment, unless you have consulted your physician and received preventative before using Silk'n Reju

Reducing the risk of injury

- Keep Silk'n Reju out of the reach of children.
- Keep Silk'n Reju out of the reach of persons with reduced physical, sensory or mental capabilities.
- Keep Silk'n Reju away from water!



Do not use or store Silk'n Reju in a wet environment (e.g. near a shower or a bathtub)!
Keep the device and its adaptor dry!

- Do not use Silk'n Reju if the ADAPTOR or the APPLICATOR is damaged!
- If the APPLICATOR is damaged, make sure it is replaced with an original part from Home Skinovations!
- The APPLICATOR contains a dedicated transformer. Never replace the APPLICATOR with another plug!

If you are unsure if Reju is safe for you to use, please consult with your doctor or dermatologist.

Possible side effects using Reju

When used according to the instructions, side effects and complications associated with use of Silk'n Reju are rare. However every cosmetic procedure, including those designed for home use, involve some degree of risk. Therefore it is important that you understand and accept the risks and complications that can occur with energy based systems designed for home use.

Silk'n Reju treatments are generally very well-tolerated. Most users sense the moderate heat applied by the device during use, with no further discomfort.

• Minor Skin Discomfort

Your skin may become slightly irritated after using Silk'n Reju. This is a typical and transient phenomenon. If the irritation does not disappear within three days, please consult a physician or a dermatologist.

• Skin Redness

Your skin may become slightly red right after using Reju. This should disappear a few hours after treatment.

Following the treatment it is recommended to apply a small quantity of the Silk'n Reju Cream.

Protect your Silk'n Reju

- Store the device in a clean area free of dust.
- Do not cover the ventilation slots during use.
- Never shake or drop Silk'n Reju or subject it to severe shocks.
- Silk'n Reju is sensitive to environmental temperature during use.
Do not expose the Reju to environmental temperatures lower than 15°C / 59°F or higher than 35°C / 95°F during use.
- Silk'n Reju can not tolerate extreme change of temperature.
When moved from a cold environment (under 10°C / 50°F) to a warm environment, avoid use for approximately 3 hours.
- Do not expose Silk'n Reju to direct sunlight or UV light for long periods of time.
- Silk'n Reju is designed for occasional home use (as opposed to exhaustive commercial or professional use). Avoid such exhaustive use.
- Avoid covering Silk'n Reju and its adaptor during recharge.

Skin Treatment the Reju Way

The benefits of Silk'n Reju

The outcome of using Silk'n Reju when following instructions correctly will be the reduction of facial fine lines, reducing pore size and improving the appearance of the face.

Best results can be expected 3-7 weeks after completing the recommended course of treatment.

Silk'n Reju facial treatment areas

Typical areas for Silk'n Reju treatment are the cheeks, forehead, neck, and around the eyes.



First Steps with Silk'n Reju

Silk'n Reju Setup

1. Remove Silk'n Reju *BASE UNIT*, *APPLICATOR*, *ADAPTOR* and other components from box.
2. Before using the Silk'n Reju for the first time, fully charge the battery. This may take up to 6 hours.

Silk'n Reju Battery Charging

- Insert the *APPLICATOR* into its cradle in the Reju *BASE UNIT*.
- Plug the *ADAPTOR PLUG* into the Base Unit socket.
- Plug the *ADAPTOR* into an electrical outlet. An orange *BATTERY – CHARGING – INDICATOR LIGHT (K)* will illuminate to indicate that your Silk'n Reju is in the process of charging.

Caution!

Never cover Silk'n Reju and its adaptor while recharging!

Note!

Silk'n Reju can not be used while recharging



- When *BATTERY – CHARGING - INDICATOR LIGHT (K)* turns off and *BATTERY – FULL - INDICATOR LIGHT (L)* illuminates, your Silk'n Reju is fully charged and ready to start.
(Recharging will be needed when *BATTERY – CHARGING - INDICATOR LIGHT (K)* illuminates again.)

Plan your Silk'n Reju treatment for best results

The efficiency of the results varies from person to person and according to the frequency of use. Some may see results sooner than others.

A single Reju session usually takes 15-20 minutes in which you should gently glide the Reju applicator slowly over the desired treatment area. During the procedure you should feel gentle heat penetrating deep into the skin.

Initial Silk'n Reju treatment course: For the first month (4 weeks), treatment should be performed 3 times a week for 15-20 minutes each time.

Silk'n Reju maintenance regime:

Use Silk'n Reju once a week for 15-20 minutes.

Treating with Silk'n Reju for the First Time

The skin should be, clean, dry and free of any powder, gel or cream.

- Remove Silk'n Reju *APPLICATOR* from its cradle in the *BASE UNIT*.
- Press the *POWER ON/OFF SWITCH*. The green *READY/TREATMENT MODE 1 - INDICATOR LIGHT (I)* will start blinking, and shortly afterwards will be lit steadily.
Your Reju is ready for treatment.
- Place the *TREATMENT SURFACE* on your skin. The *LED* light will automatically activate.
- Apply light pressure while slowly moving the Reju *APPLICATOR* in small circular motions. The *TREATMENT SURFACE* temperature will slowly increase.
- Perform the treatment for 5-7 minutes on each area.
- It is recommended to apply a small amount of a moisturizing cream on the treated skin areas after each session and to massage it gently in circular motions.



What to expect when treating with Silk'n Reju?

During a Silk'n Reju session it is normal to experience and feel:

- **A Sensation of Warmth** – During the Silk'n Reju session you may feel slight warmth. This is normal and should disappear up to an hour after the session.
- **Some Mild Red Color** – After each Silk'n Reju session you may notice a slight reddening of the skin (mild erythema) over the treated area. This is normal and should fade within an hour or two.

Best results can be expected 3-7 weeks after completing the recommended course of treatment.

After treating with Silk'n Reju

- When Reju session has been completed turn Reju off by pressing the *POWER ON/OFF SWITCH*.
- After each session it is recommended that you clean your Reju device, especially the *APPLICATOR TIP* (See Section: "Cleaning Reju").
- Place the *APPLICATOR* in its *BASE UNIT* cradle to recharge the batteries (See Section: "Silk'n Reju Battery Charging").

Maintenance of Silk'n Reju

Cleaning the Silk'n Reju

After each session, it is recommended to clean your Silk'n Reju device:

- Turn Silk'n Reju off by pressing the *POWER ON/OFF SWITCH*.
- Unplug Silk'n Reju from the electrical outlet
- Allow cooling down of its *TREATMENT SURFACE*.



- Use a dry, clean cloth and a specially formulated cleaner for electronic equipment to gently wipe the device and especially its *TREATMENT SURFACE*.



- Never immerse Silk'n Reju or any of its parts in water!

- Never clean the Reju or any of its parts under the tap or in the dishwasher.

Storage

- Store the clean device and all its parts in a dry place, at a temperature of between -10°C to 60°C / 14°F to 140°F

To protect the environment

At the end of its life, do not discard Silk'n Reju with the normal household waste. Dispose of it at an official recycling collection point.

Remove batteries and dispose at an official collection point for batteries. Please call your local recycling depot for locations and proper recycling instructions.

Removing the rechargeable batteries

Remove the rechargeable batteries only if they are completely empty and cannot be recharged.



Undo the screw at the bottom panel of the *APPLICATOR*.



Lift and remove the top panel of the *APPLICATOR*.



Push the three snaps and pull out the electronics card.



Cut the tow wires that connect the rechargeable batteries pack to the *APPLICATOR*.



Cut the tow wires that connect the electronics card to the lower panel of the *APPLICATOR*.

Troubleshooting

This chapter summarizes the most common problems you could encounter. If you are unable to solve the problem with the information below, contact your local Silk'n Customer Service Center for further support.

Problem	Possible Cause	Solution
Silk'n Reju does not work	The rechargeable batteries are empty	Charge the <i>APPLICATOR</i>
	The <i>APPLICATOR</i> is broken	Contact your local Silk'n customer Service center
Silk'n Reju suddenly switches off	The rechargeable batteries are empty	Charge the <i>APPLICATOR</i>
Device temperature rises during recharge	This is normal	No action is required
BATTERY – CHARGING - INDICATOR LIGHT (K) does not turn on when connecting the PLUG into the BASE UNIT SOCKET	The <i>ADAPTOR</i> is not plugged in the electrical outlet	Insert the <i>PLUG</i> into the <i>BASE UNIT SOCKET</i>
		Connect another appliance to the wall socket to check if the electrical outlet is live
		If the <i>APPLICATOR</i> still does not charge, contact your local Silk'n Customer Service Center
	The <i>ADAPTOR</i> is not inserted properly into the electrical outlet and/or the <i>PLUG</i> is not inserted properly into the <i>BASE UNIT SOCKET</i>	Make sure the <i>APPLICATOR</i> is inserted properly into the electrical outlet and/ or <i>PLUG</i> is inserted properly into the <i>BASE UNIT SOCKET</i>
	The <i>APPLICATOR</i> is broken	Contact your local Silk'n Customer Service Center
APPLICATOR temperature rises during use	This is normal	No action is required
READY/TREATMENT MODE - INDICATOR LIGHT (I) does not illuminate when pressing the ON/OFF SWITCH	The <i>APPLICATOR</i> needs to be reset	Press the <i>ON/OFF SWITCH</i> twice to switch the <i>Applicator</i> off and on again. If <i>READY/ TREATMENT MODE - INDICATOR LIGHT (I)</i> does not illuminate again, contact your local Silk'n Customer Service Center

The Yellow SYSTEM ERROR - INDICATOR LIGHT (j) is lit.	Malfunction of the APPLICATOR	Press the ON/OFF SWITCH twice to switch the Applicator off and on again. If SYSTEM ERROR - INDICATOR LIGHT (j) illuminates again, contact your local Silk'n customer Service center
The Yellow SYSTEM ERROR - INDICATOR LIGHT (j) is blinking	The rechargeable batteries are empty	Charge the ADAPTOR
The Orange BATTERY – CHARGING and Green BATTERY – FULL - INDICATOR LIGHTS are illuminated alternately	Malfunction of the ADAPTOR	Disconnect the ADAPTOR from the electrical outlet, press the ON/OFF SWITCH twice to switch the APPLICATOR on and off again. Connect the Adaptor to the wall socket, If the problem persists, contact your local Silk'n Customer Service Center
	Malfunction of the APPLICATOR	Press the ON/OFF SWITCH twice to switch the Applicator off and on again. If the problem persists, contact your local Silk'n Customer Service Center

Customer Service

For USA and Canada customers:

For more information about Silk'n™ products please visit www.silkn.com.

If your device is broken, damaged, in need of repair; or for any other user assistance, please contact Silk'n™ Customer Service: **1-877-DO-SILKN / 1-877-367-4556**, contact@silkn.com

International Customers:


For more information about Silk'n products or to find a list of International Silk'n websites please visit www.silkn.eu or www.silkn.com.

If your device is broken, damaged, in need of repair; or for any other user assistance, please contact your local Silk'n™ Customer Service or write to: info@silkn.com.


Manufacturer:

Home Skinovations Ltd.
 Apolo Building, Shaar Yokneam POB 533
 Yokneam 20692
 ISRAEL





Specifications

Rated voltage (Adaptor)	100V - 240V
	Model: UE08WCP-105040SPA
Rated frequency (Adaptor)	50Hz – 60Hz
Rated current (Adaptor)	400 mA
Applied part	 Type BF
Protecting rating	IP20 (EN60529)
Storage condition	Temperature: -10°C to 60°C / 14°F to 140°F
	Relative humidity: 5% to 95%
Operating condition	Temperature: 15°C to 35°C / 59°F to 95°F
	Relative humidity: 10% to 80%
	Atmospheric pressure: 700 to 1060 hPa
Light source	Light Emitting Diodes
	Wavelengths peak emission wavelength: 645nm
	Intensity power: 250mW ±35%
Lithium ion battery	2 × 3.7V, 750±50mAh
Max. charging time	180 minutes


Labeling

HomeSkinovations  Appolo Building
Shaar Yokneam, P.O.Box 533
20692 Yokneam, Israel


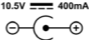
MODEL : Reju
PART NO. : AS1018321
MANUFACTURED : 2610





For use only with power adaptor specified in user's manual
Pour l'usage seulement avec l'adaptateur de puissance indiqué en manuel d'utilisateur

000100001


Made in ISRAEL


 10.5V  400mA

LB1019521

Degree of protection against electric shock type BF applied part	
WEEE - Waste Electrical and Electronic Equipment	
CSA Mark	
CE Mark	

Degree of protection against ingress water: ordinary

This device is not suitable for use in the presence of flammable anesthetic mixture with air or with Oxygen or Nitrous Oxide.

The symbol  on the label affixed to this device means "Attention, consult accompanying documents".

EC Declaration of Conformity

We Home Skinovations Ltd. declare under our sole responsibility that the home Skin treatment products

Brand Name Home Skinovations

Model name Model Number

Silk'n Reju AS1018321

Are fully in conformity with the essential requirements of EMC Directive 2004/108/EC and Low Voltage Directive 2006/95/EC.

This declaration is based on the full compliance of the products with the following European standards:

EMC: EN 55022:2006
EN 55024:1994;(A1:2001)(A2:2003)
EN 61000-3-2:2006; (A2:2005),
EN 61000-3-3:1995; (A1:2001)

Electrical Safety: EN 60335-1:(2002),
EN 60335-2-23:(19963)

Authorized Representative: Duotec Medical GmbH. Dr. Kurt Huber Str. 6, D-82031 Grünwald, Germany
Contact person: Hans Edel

Manufacturer: Home Skinovations Ltd. 4 Hacharoshet St., Kfar Saba P.O.B. 7028, Israel, 44640, Tel: (972) 9 7678807 | Fax: (972) 9 7678808



Gabi Lavi
General Manager
Yokneam, Israel
July 12, 2009



Hans Edel
Europe
July 12, 2009


HomeSkinovations

Apolo Building, ShaarYokneam POB 533, Yokneam 20692, ISRAEL www.silkkn.com
tel. +972-4-9097470, fax +972-4-9097471, info@silkkn.com, contact@silkkn.com



Apolo Building, Shaar Yokneam POB 533
Yokneam 20692, ISRAEL
tel. +972-4-9097470, fax +972-4-9097471
info@silkn.com, contact@silkn.com
www.silkn.com



PB101950A